




# Cap-Shure® Cotton Tipped Applicator

Polystyrene handle

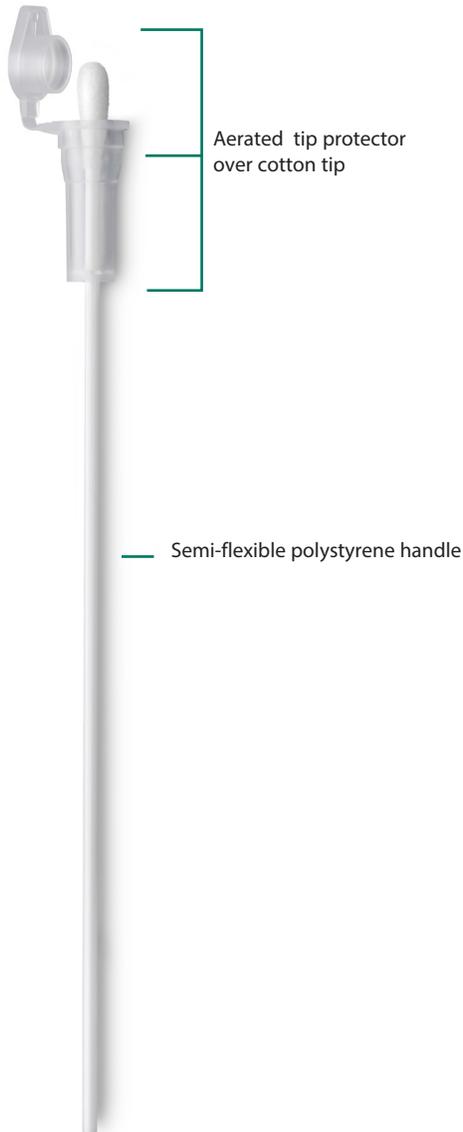

Aerated tip protector  
over cotton tip

Semi-flexible polystyrene handle

\* Not actual size

## Description

This absorbent cotton tipped applicator is produced with high quality USP grade cotton which is carded and coiled onsite to maintain consistency. It is bonded securely to a polystyrene handle by an aqueous based adhesive. Aerated tip protector designed for convenient, reliable specimen, DNA or evidence collection. Comes as a unit with a reclosable, ventilated cap for air drying and transporting samples.

## Dimensions

6.125" Overall Length (155.57 mm)

### Tip

- 0.203" Dia. (5.156 mm)
- 0.688" Length (17.47 mm)

### Handle

- 0.100" Dia. (2.54 mm)
- 5.875" Length (149.225 mm)

## Qualities

- Medical grade quality
- Tip suspended away from sides
- Can lay flat to air dry
- Individually wrapped
- Made in USA

## Applications

Diagnostic /Forensics & DNA

- Specimen collection
- DNA collection
- Evidence collection
- Paternity testing
- For use in the field or home testing
- Buccal swabbing

Ref #	Packaging
25-806 1PC EC Sterile	10/50/1 (10 boxes per case, 50 individually wrapped packages per box, 1 p/pkg)



**Puritan®**